


Likorall™ Overhead Lifts – a name recognized in hospitals, nursing homes, rehab facilities and home care. A product known for ease of use, well-engineered design and a genuine regard for safety, has a new member of the family: the Likorall™ 200 Overhead Lift.

The Likorall™ 200 lift uses the same proven features and technology designed for safety as other overhead lifts in the Likorall™ series but with a simplified design. This new design makes the Likorall™ 200 lift an uncomplicated and affordable lift motor to solve your daily home lifting needs.

Key features at a glance

- Ideal for basic straight rail or straight rail with a curve design, or for use on Liko® FreeSpan™ or FreeStand™ products
- Easy to use integrated two button hand control, with spring clamp for quick access
- QLink included for use with Liko's quick release accessories
- Hill-Rom's patented safety drum
- Electrical emergency lowering with the pull of a cord

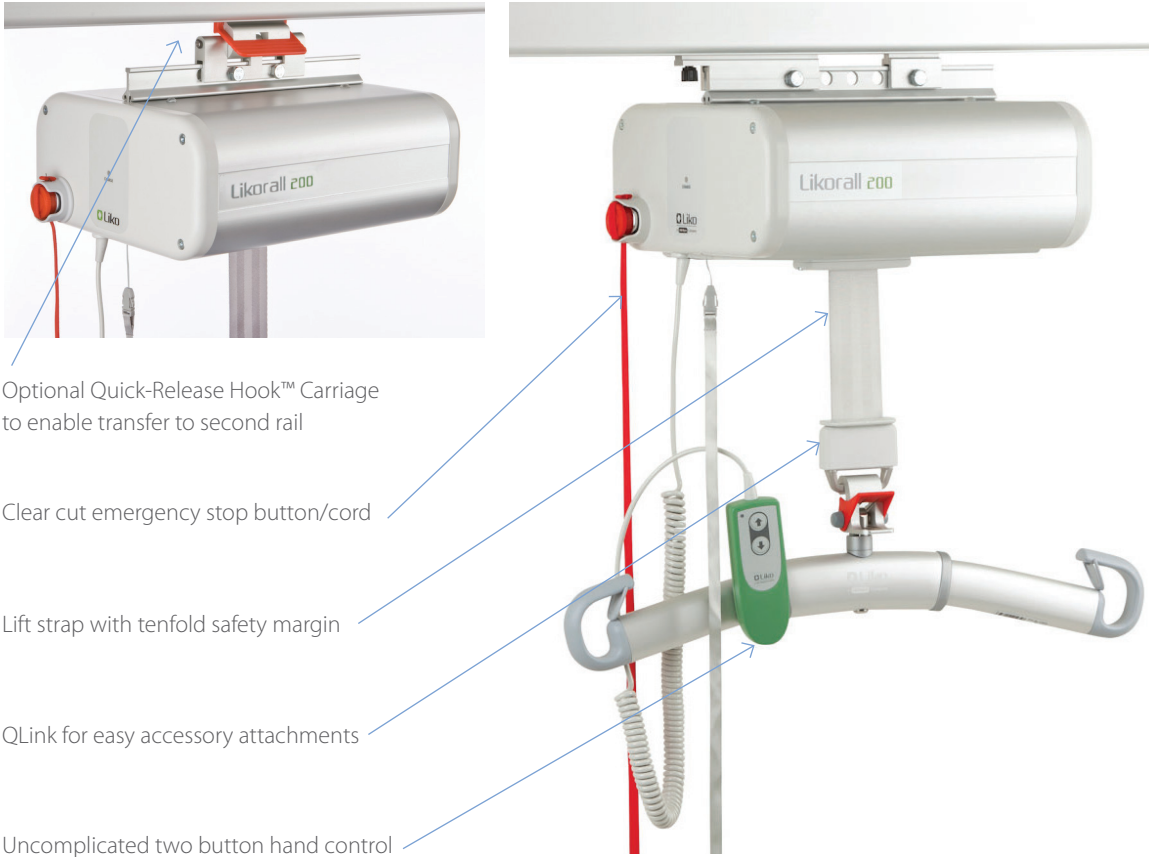

Optional Quick-Release Hook™ Carriage
to enable transfer to second rail

Clear cut emergency stop button/cord

Lift strap with tenfold safety margin

QLink for easy accessory attachments

Uncomplicated two button hand control

Liko® Likorall™ 200 Overhead Lift

Max Lifting Capacity	440 lbs
Lifting Speed	2 inch/s
Lifting interval	80.7 inches
Lift Unit Weight	27.5 lbs
Lift Unit Dimensions	13.4 x 9.8 x 6.5 inch (L x W x H)
Hand Control Cable length	34 inches
Electrical Data	24 V, 12 A

The Likorall™ lift is equipped with a Single Fault Safety (SFS) drum. This patented design provides protection against uncontrolled lowering.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2014 Liko R&D AB. ALL RIGHTS RESERVED.

ORDER NUMBER 186110 rev 1 25-FEB-2014 ENG – US

Zinc Number US-PA-0114-0036

www.hill-rom.com

USA 888-545-6671